

Image Quality.
Ultimate Access.
Versatility.

Stryker Prime X™ X-ray Stretcher

Stryker stretchers have been the industry benchmark for decades, and Stryker's Prime Series Stretchers continues the tradition.

The Stryker Prime Solution, Now with X-Ray Capability

Advanced mobility, safety and efficiency combine with the Clearview Technology platform.

The Prime X X-Ray Stretcher is the next generation of Prime Series Stretchers, loaded with all of the features you already love — now, with X-ray capability. Stryker's exclusive Clearview Technology platform provides similar image quality to that of a stationary X-ray table, on the Prime Series Stretcher mobile platform.

- 360 degree access for cassette placement
- 52% larger imaging window
- 700 lb. weight capacity
- Integrated scale
- Power-washable design

The results? Increased safety; incredible efficiency.

Pneumatic Backrest

Allows caregivers to raise the backrest without the risk of damaging expensive cassettes and detectors.

Clearview Technology™ Platform

Structural foam encased by two layers of ABS composite offers comparable image quality to X-ray tables.

Easy ID

The Prime X X-ray Stretcher exclusively utilizes a red color scheme on its handles, side-rails, bumpers and accent rings so that your X-ray stretcher is easily identified.

Electric Lift

With electric lift, simply touch the pedal to raise the litter. Available on Fifth Wheel and Big Wheel Advanced Mobility models.



Advanced Mobility

Choose from retractable Fifth Wheel traction and cornering, Big Wheel can reduce start-up forces by up to 60%, or Zoom® Motorized eliminates manual pushing.

Prime X X-ray Stretcher Accessories



Lateral Cassette Holder




Upright Cassette Holder



Foot End Push Handles
(1105x only)

The Prime X X-ray Stretcher offers the safety and efficiency of the Prime Series Stretcher, now with the Clearview Technology platform.



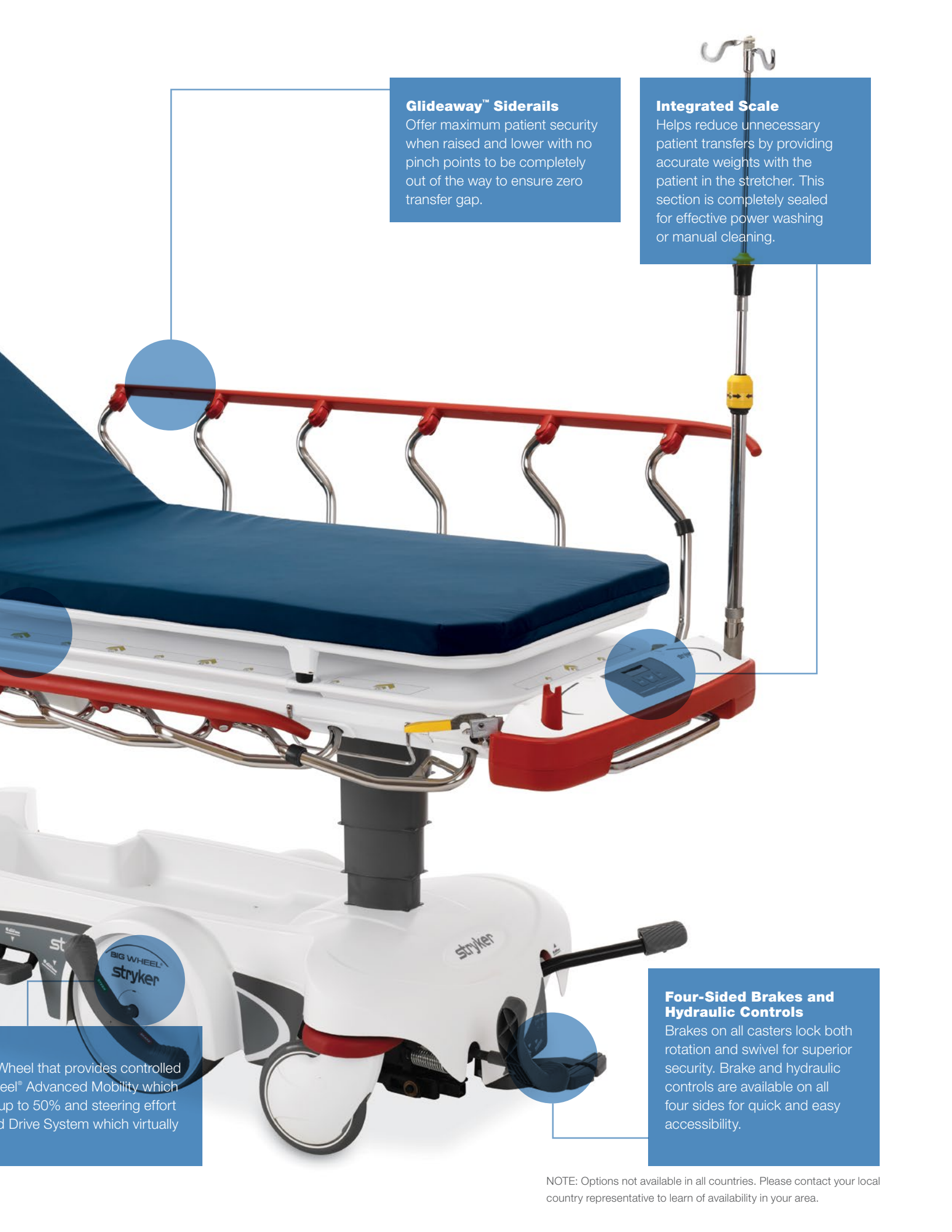

Glideaway™ Siderails

Offer maximum patient security when raised and lower with no pinch points to be completely out of the way to ensure zero transfer gap.

Integrated Scale

Helps reduce unnecessary patient transfers by providing accurate weights with the patient in the stretcher. This section is completely sealed for effective power washing or manual cleaning.

Four-Sided Brakes and Hydraulic Controls

Brakes on all casters lock both rotation and swivel for superior security. Brake and hydraulic controls are available on all four sides for quick and easy accessibility.

Wheel that provides controlled
heel® Advanced Mobility which
up to 50% and steering effort
d Drive System which virtually

NOTE: Options not available in all countries. Please contact your local country representative to learn of availability in your area.

The Prime X X-ray Stretcher allows you to **adjust the cassette to the patient** rather than the patient to the cassette.

The unique dual-deck design allows positioning of cassettes at any point under the patient, from head to foot and side-to-side. Open architecture accommodates virtually any size of cassette.



Easy to Raise Foot Section offers enhanced cleanability through convenient access to the cassette platform.



Open Cassette Platform and a full 360° of access allows positioning of cassettes from head end to foot end.



Beveled Edges and Visual Alignment Guides make positioning the cassette to the patient simple.



Mobile X-ray platform helps reduce patient transfers, offering a versatile solution throughout the hospital.

Reconstructive

Hips
Knees
Trauma & Extremities
Joint Preservation
Orthobiologics

Medical & Surgical

Power Tools & Surgical Accessories
Image Guided Navigation
Endoscopy & Arthroscopy
Integrated Communications
Beds, Stretchers & EMS
Sustainability Solutions

Neurotechnology & Spine

Craniomaxillofacial
Interventional Spine
Neurosurgical, Spine & ENT
Neurovascular
Spinal Implants

3800 E. Centre Ave.
Portage, MI 49002 U.S.A.
t: 269 329 2100 f: 269 329 2311
toll free: 800 327 0770

www.stryker.com

The information presented in this brochure is intended to demonstrate Stryker products. Always refer to the package insert, product label and/or user instructions before using any Stryker product. Products may not be available in all markets. Product availability is subject to the regulatory or medical practices that govern individual markets. Please contact your Stryker Account Manager if you have questions about the availability of Stryker products in your area.

Products referenced with [™] designation are trademarks of Stryker.
Products referenced with [®] designation are registered trademarks of Stryker.

Mkt Lit-625 13 DEC 2011 Rev B.0
MHD/PS **3M** 1714, SPH255

Copyright © 2014 Stryker
Printed in U.S.A.