

# SofCare® Stretcher Overlay


Convenient. Complete.  
Effective.



The SofCare Stretcher Overlay allows hospitals to provide a **continuum of preventative skin** care from the moment a patient is admitted to the time the patient is discharged. The overlay is composed of three distinct layers containing **hundreds of compensating air cells** that are interconnected through a series of air channels. As the patient lays on the overlay, the cells exchange air and adjust for support. The patient's weight is redistributed over the entire surface, mimicking the buoyancy of fluid flotation. After use, the overlay is easily disposed of.



**Advanced pressure redistribution overlay is designed specifically for stretchers for preventative skin care or to begin treatment of pressure ulcers immediately.**

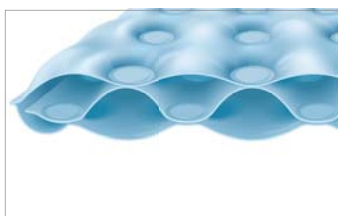

**Three layer technology with unique pocket and button design allows for air adjustment and support while helping ensure that patients don't bottom out.**



**Multiple pump options are available for convenience. Choose from the hand pump, Electric Inflator Pump, or the Companion Pump, which keeps continuous air in the overlay.**

**Single patient use** promotes infection control

**Pocket and button design** keeps patients from bottoming out

**Continuous air adjustment** provides patient with comfort and support

**Multiple pump options** allow the hospital to choose what works best for their needs

**CPR Release** is a quick and convenient pull in case of emergency

Available in **26-inch and 30-inch** models

# SofCare<sup>®</sup> Stretcher Overlay



**Medical**

3800 E. Centre Ave.  
 Portage, MI 49002 USA  
 t: 269 329 2100 f: 269 329 2311  
 toll free: 800 787 9537

[www.stryker.com](http://www.stryker.com)

**Standard Features**

- Adjustable head and side straps
- Hand welds
- Inflation valve
- CPR Release

**Optional Features**

- Manual Hand Pump
- Electric Inflator Pump
- Companion Pump

**Specifications**

Model Number	SC2775	SC3075
Overall Width	27" (70 cm)	30" (76 cm)
Overall Length	75" (191 cm)	75" (191 cm)
Weight Capacity	500 lb (226.8 kg)	500 lb (226.8 kg)

**Warranty**

- Manufacturing defect only

*Stryker reserves the right to change specifications without notice.*